


BAT BRIGHTNESS ACUITY TESTER

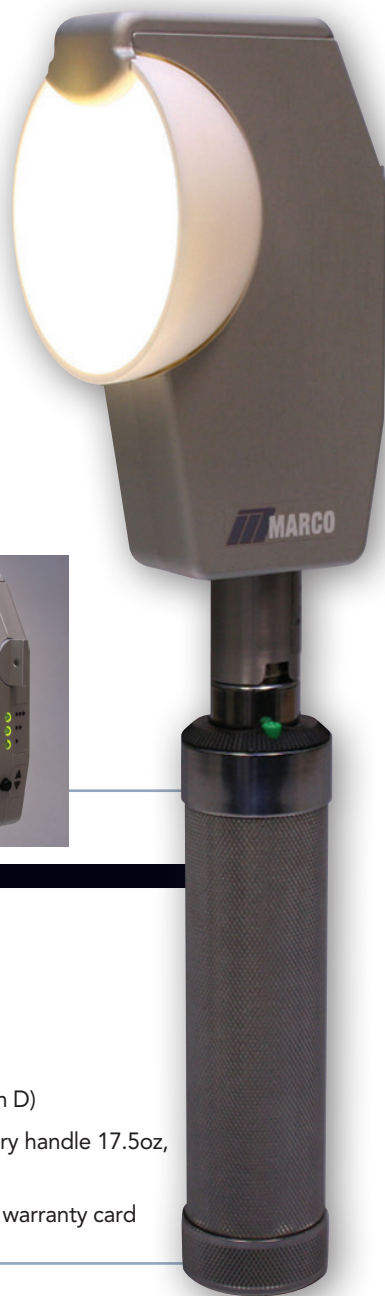

BAT Brightness Acuity Tester

FEATURES/SPECIFICATIONS

Light output:	HIGH = 400 foot-lamberts MEDIUM = 100 foot-lamberts LOW = 12 foot-lamberts
Power Input:	3.5VDC 750ma
Size:	5" H X 2.5" W X 2.5" D (128mm H X 64mm W X 64mm D)
Weight:	0.5 lb. (.23 kg) (With WA battery handle 17.5oz, 1lb 1.5oz, 496g)
Standard Accessories:	BAT head, instruction manual, warranty card

- Recognized industry-standard glare testing device effectively measures glare without leaving the exam lane. Can measure all patients with any standard acuity chart.
- NEW modular, solid-state design accommodates most standard (twist-lock) rechargeable instrument handles. BAT recharging is easily accomplished via an instrument stand well, or standard desk-top charger.
- NEW state-of-the-art LED illumination system provides a better, more efficient light source.
- Three variable light intensity settings simulate various glare conditions.
- Evaluates tinted lenses, sunglasses and anti-reflectives to quantify improvement in glare disorders.
- Effectively evaluates Macular photo stress testing.
- Simple, very easy to use.

Protect your warranty. Buy your Marco product only from an authorized Marco distributor. Call Marco for the name of the distributor nearest you.

