

MARCO ENCORE AUTOMATIC CHAIR





MARCO Encore Automatic Chair

With the touch of a button, the Marco Encore Automatic Chair quickly reclines from an upright position to a fully horizontal, flat-table mode—or anywhere in between. Another button instantly returns the chair to its upright position. Marco's exclusive hydraulic system ensures silky smooth elevation and recline. This unique, fail-safe hydraulic mechanism also enables the chair to descend to an extremely low position, yet the high-end elevation remains uncompromised. Backlit flat overlay switches on both sides of the chair are easy to see in darkened rooms. The Marco Encore Automatic Chair is the ideal choice for those who prefer instant, automatic recline functions. Electrical requirements: 110V AC, 50/60 Hz, 7 amps maximum.

- Easy-to-use backlit fingertip switches on both sides control chair elevation, recline and return functions. Programmability allows the touch of a switch to instantly recline the chair to form a bed or any position desired.
- Recline, as well as elevation, is controlled by Marco's exclusive field-proven hydraulic system for smooth, quiet adjustments.
- Corded foot-pedal allows auxiliary control of elevation, recline and automatic return.
- Can be interfaced with Marco's Deluxe or Bravo Stands to form a wheelchair-accessible system.
- Standard single-handed headrest with unique lock indicator assures the headrest is firmly and securely in position. Optional pillow headrest and surgical package are also available for minor surgical procedures.
- Scuff-resistant, high-impact covers on arms and chair base are easily cleaned.
- Standard Marco ten-step, hand-finished painting process.




Single-handed headrest



Easy, "dial-in" programming control with safety switch to disable chair



The Encore Automatic can be instantly transformed into a surgical/minor procedure chair.

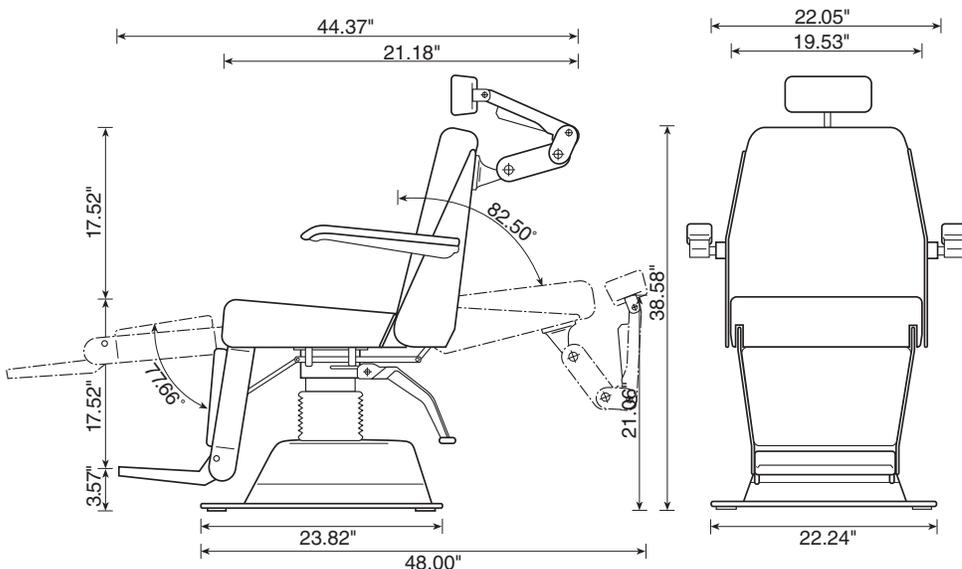

Backlit flat overlay switches are easy to see in a darkened room.



800.874.5274
www.marco.com