



CHAIRS & STANDS



MARCO CHAIRS & STANDS

The Difference is Marco

MARCO E-Z TILT

The Marco E-Z Tilt Chair combines the most simplified manual reclining motion with a new modern, sophisticated appearance. The Marco E-Z Tilt also incorporates the only double hydraulic cylinder system in the industry, creating an extended range of elevation, from extra low to extra high, helping to accommodate a full range of patient and practitioner sizes. Convenient controls for chair reclining, elevation, rotation, and headrest adjustments are fingertip and foot operated.



EXCLUSIVE TELESCOPING
HYDRAULIC CYLINDER SYSTEM



MARCO ENCORE AUTOMATIC


The Marco Encore Automatic Chair features an exclusive hydraulic system that controls both the elevation and recline motion. Easy-to-use backlit fingertip switches control all functions of the chair, as well as programmable recline positions, auto-recline/return, and safety functions. A corded foot pedal and single-handed operation headrest are also standard. Optional surgical accessories are available for minor surgical procedures.



EASY, "DIAL-IN"
PROGRAMMING CONTROL
WITH SAFETY SWITCH TO
DISABLE CHAIR OPERATION

MARCO ENCORE MANUAL

Considered our standard, everyday exam chair, the Marco Encore Manual Chair incorporates a simple, pneumatic-assisted recline mechanism for easy patient adjustments. Convenient fingertip switches on the side of the chair instantly control chair elevation and auto return. A corded foot pedal and single-handed operation headrest are also standard.

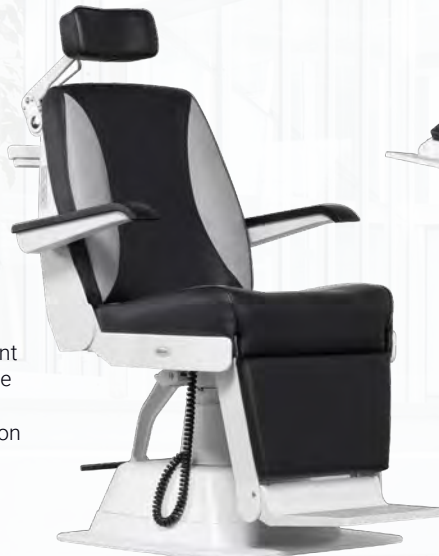



MARCO BRAVO 2

The Marco Bravo 2 Exam Chair offers the same great value of the original plus enhanced functionality and features. In addition to a smooth, quiet hydraulic lift, an effortless tilting recline mechanism is easy to access and keeps the patient in a comfortable position. With its modern, ergonomic design, and robust engineering, the Bravo 2 chair is made to look good and last.




BRAVO 2, EFFORTLESS TILT



MARCO DELUXE STAND

A standard in the industry, the modern Marco Deluxe Instrument Stand meets the highest demands of function, durability, and cosmetic appeal. Innovative features include a re-designed lower slit lamp arm with multiple integrated adjustment functions, a modified refractor arm and instrument pole that provides convenient routing of cables, and a simplified membrane console panel.

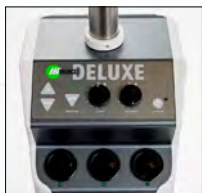

DELUXE MEMBRANE CONSOLE PANEL



MARCO BRAVO 2 STAND | NEW

The Bravo 2 Instrument Stand from Marco combines a fresh, new updated design with our high standards of performance, function, and durability. With a strong emphasis on simplicity and economy, the Bravo 2 stand provides the perfect blend of everyday efficiency and value.



BRAVO 2 MEMBRANE CONSOLE PANEL



TECHNICAL SPECIFICATIONS

EZ TILT EXAMINATION CHAIR

- » Base: 22.25" x 23.85"
- » Lift Capacity: Approx 400 lbs
- » Ground to Seat Height: 20.5" low – 33" high
- » Rotation: 270° degrees
- » Recline Mechanism: Tilting recline
- » Input Voltage: 120VAC 60HZ
- » Fuses: 8A 250V
- » Available Colors: Black, Charcoal, Navy, Silver, Silver & Black two tone, Silver & Navy two tone

BRAVO 2 EXAMINATION CHAIR

- » Base: 22" x 23.25"
- » Lift Capacity: Approx. 350 lbs
- » Ground to Seat Height: 21.5" low – 30" high
- » Rotation: 270° degrees
- » Recline Mechanism: Tilting recline
- » Input Voltage: 100–120VAC 60HZ
- » Fuses: 8A
- » Available Colors: Black, Charcoal, Navy

ENCORE MANUAL EXAMINATION CHAIR

- » Base: 22.25" x 23.85"
- » Lift Capacity: Approx. 375 lbs
- » Ground to Seat Height: 21" low – 32" high
- » Rotation: 270° degrees
- » Recline Mechanism: Manual Lay Back
- » Input Voltage: 120VAC 60HZ
- » Fuses: 8A 250V
- » Available Colors: Black, Charcoal, Navy, Silver, Silver & Black two tone, Silver & Navy two tone

ENCORE AUTOMATIC EXAMINATION CHAIR

- » Base: 22.25" x 23.85"
- » Lift Capacity: Approx. 375 lbs
- » Ground to Seat Height: 21" low – 32" high
- » Rotation: 270° degrees
- » Recline Mechanism: Auto Lay Back
- » Input Voltage: 120VAC 60HZ
- » Fuses: 8A 250V
- » Available Colors: Black, Charcoal, Navy, Silver, Silver & Black two tone, Silver & Navy two tone

DELUXE INSTRUMENT STAND

UNIT

- » Base: 21.25" x 21.25"
- » Console Height: 33.5"
- » Power: AC 120V (50/60 Hz)

SLIT LAMP ARM

- » Capacity: Up to 60 lbs
- » Range: 13"

REFRACTOR ARM

- » Capacity: Up to 12 lbs
- » Range: 8"

BRAVO 2 INSTRUMENT STAND

UNIT

- » Base: 22" x 23.25"
- » Height Console: 33.5"
- » Height Pole: 68"

SLIT LAMP ARM

- » Capacity: Up to 60 lbs
- » Range: 13.5"

REFRACTOR ARM

- » Capacity: Up to 20 lbs
- » Range: 12.5"

Contact Marco First. MARCO.COM | 800.874.5274

Marco Healthcare a division of Advancing Eyecare™, 11825 Central Parkway, Jacksonville, Florida 32224